

For your patients!



TOPI-CLICK[®]

DISPENSING INSTRUCTIONS*

Read all of this information and instructions carefully before dispensing from the Topi-Click[™].

1. Pop off the cover by placing thumb in groove on the front where it says REFILL and press upward on the cover's thumb tab.
2. Turn the base slowly as shown below in Fig. 2 to dispense the specified ____ # of clicks per topical provider's instructions. You may refer to the bottom for reminder instructions
3. Apply by rubbing applicator pad to area of skin as directed by topical provider.
4. Replace cover firmly and tightly.
5. When the bottom line (see Fig. 5 below) reaches the engraved line above the word "Refill," it may be time to call in for any refills because there are approximately 12 clicks left.

***The Refill Indicator:** There are approximately (depending on viscosities) the following numbers of clicks left when the very bottom of the topical elevator reaches the line engravings on the front; refer to Fig. 5 below.

➤ looking at the top front there are these engravings:

____ = 12 clicks left (first line from the top)


Refill

____ = 32 clicks left (second line from the top)



The Topi-Click[®] is a one-time use container and not refillable.

Downloadable at www.Topi-Click.com

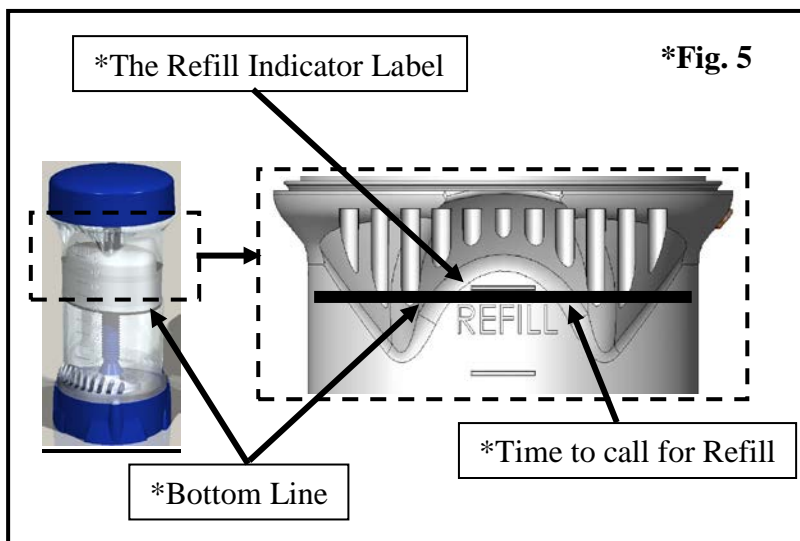

The Topi-Click[®] product may be improved over time based on customer feedback on usage. Custom Rx TDA, LLC and/or the patent holders are not liable for any damages arising from the products dispensed in the Topi-Click[®] or its improper use.



Protected by U.S. Patent No. 7,213,994.

TOPI-CLICK

1309 Crooked Branch Trail
Woodstock, GA 30189-6830

www.Topi-Click.com
770-827-0472



All items © 2006-2011 by Custom Rx TDA, LLC.
All rights reserved Ver: 72010